

USER'S MANUAL
GLASS PERSONAL SCALE

Code: SCALE-PRS-BLOOM-WL



Thank you for choosing **WELL**. Please read carefully the following instructions and keep them within reach.

Specification

- Capacity: 180kg
- Division: 100g
- Auto zero & auto off function
- Low battery and overload indication
- Power: 1x 3V (CR2032 button cell) included

CAUTION


1. Keep the scale off water



2. Keep the scale with care



3. Do not put it on uneven floor while weighing

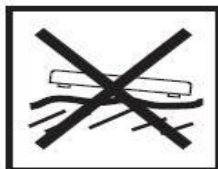

4. Do not step on the edge of the platform



5. Do not step on the scale with wet feet

6. Do not use the scale under electromagnetic interference condition

7. Please take the batteries out if the scale is not used for a long time

8. Do not attempt to open the scale or to remove any components. Servicing should be done by qualified technicians only.

BEFORE USING SCALE


1. Remove all packaging materials.

2. Open the battery cover, insert 1XCR2032 batteries to the battery case. Then close the battery cover. Or just remove insulating strip.

USER'S MANUAL
GLASS PERSONAL SCALE

Code: SCALE-PRS-BLOOM-WL

Well®

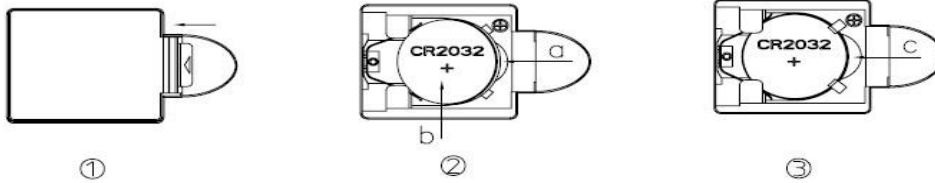


3. Set scale on floor, Surface must be level for accurate operation.
4. Push C position to take out batteries easily.


To Measure Weight

Before your first using, please follow below step:

1. Place the scale on a flat horizontal surface, and the place must be hard .This condition is obligatory – it allows to receive an exact results . Then directly stand on the scale with your foot to turn it on..The display lights up. As shown in the follow picture

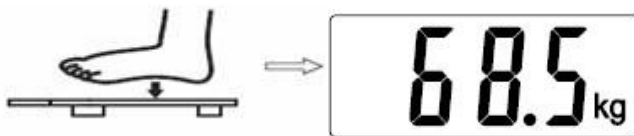

2. The factory default unit is “kg” , if not applicable, please set switch on the back of scale to choose the needed unit lb or st.
3. Auto onto the scale

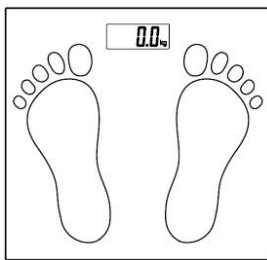

4. The LCD displays your Weight: As shown in the follow picture



The scale will turn off automatically after 10 -15 seconds without any operations.

USER'S MANUAL
GLASS PERSONAL SCALE

Code: SCALE-PRS-BLOOM-WL

Well

NOTE


1. "C" will appear if the scale is in unstable conditions, it means it's not the real weight, you need to be weighted again.



2. "ERR" will appear if the scale is overloaded (Over 180kg/400lb).



3. "Lo" will appear if battery is low and needs to be replaced.

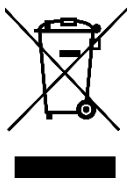

Waste electrical and electronic equipment are a special waste category, collection, storage, transport, treatment and recycling are important because they can avoid environmental pollution and are harmful to health. Submitting waste electrical and electronic equipment to special collection centers makes the waste to be recycled properly and protecting the environment. Do not forget! Each electric appliance that arrives at the landfill, the field, pollutes the environment!



Importer & distributor:
SC VITACOM ELECTRONICS SRL

CIF: RO 214527

Tel. 0264-438401*

sales@vitacom.ro, www.vitacom.ro